

VEYA™

EFFORTLESS AUTOMATION.

Clinical Suite.

Take the first steps toward clinical workflow automation with Veya, built on Tecan's trusted modular platform, proven in clinical markets.


Pipetting with precision.

Four-channel, air-restricted pipetting enables precise and consistent multi-sample processing for reliable and reproducible results.

Regulatory compliance.

GX assurance software ensures data integrity and comprehensive audit trails, including support for electronic signature.

Worktable.

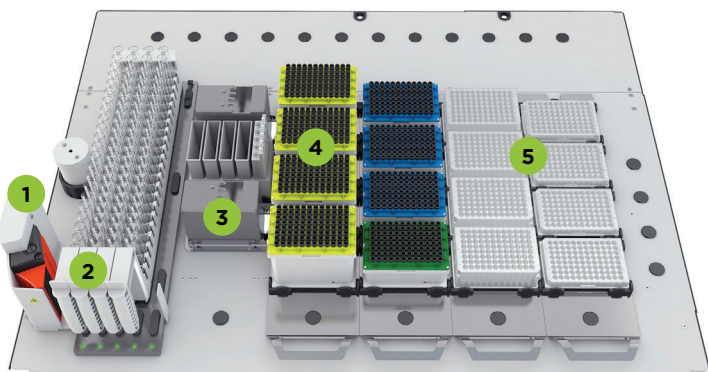

Example worktable for the Veya Clinical Suite.

To learn more, please visit www.tecan.com/veya

*Class 1 laser product pursuant to IEC 60825-1:2014

**EU Regulation 2017/746

The type of specimen and specific diagnostic protocol is defined and validated by the user for their specified clinical assay. Contact the Tecan sales organization for the local regulatory status of the product.



Flexible capacity.

Reformat up to 120 samples from tubes to 96-well plate format.

Sample traceability.

Integrated loading ID and automatic barcode scanning guarantee sample traceability, eliminating manual errors and saving time.

Ready for regulated environments.

Compliant to IVDR class A** and USA class 1 IVD, 510(k) exempt.

- 1 Loading ID barcode scanner*
- 2 Runner segments
- 3 Integrated waste management system (tip-drop station)
- 4 Disposable filtered tips in ANSI/SLAS boxes (SBS Combi format)
- 5 Deck segments for microplates
- + Includes independent flexible four-channel arm with Air Restricted Pipetting (ARP) technology (shown in full system view)

LEARN MORE



403128 V1.0, 2025-01

